

ROOT TEASERS and ELEVATORS

Heidbrink Root Tip Teasers

- Unique instruments for teasing out root tips from alveoli and canals.
- Heidbrinks have a special curved point for insertion into the ligament space for gentle lifting.
- Each instrument has a right and left side to work mesially and distally.
- Heidbrinks are NOT ELEVATORS - the fine tips will fracture; use with a gentle lifting and pulling force.
- Instruments with bent or dull tips should be replaced.

Cat. No. 875-001 Heidbrink #1 - 9.50 each
 Cat. No. 875-002 Heidbrink #2 - 9.50 each
 Cat. No. 875-003 Heidbrink #3 - 9.50 each

SPECIAL

1 set 3 Heidbrink
Cat. No. 875-000
One of each size
 List \$28.50 for \$25

SIZE	WIDTH	CAT. NO.
1	1.2 mm	875-001
2	1.5 mm	875-002
3	1.8 mm	875-003




Kopp Elevators

- Ultra thin tip to enter ligament space.
- Razor-sharp tip also used to pierce root tip for extraction.
- Each instrument has a right and left side to work both mesially and distally.
- Place elevator in ligament space and use a circular motion to extract root tip.
- Bent or dull tip instruments should be replaced.

Cat. No. 875-011 Kopp Left - 21.00 each
 Cat. No. 875-012 Kopp Right - 21.00 each
 Cat. No. 875-013 Kopp Left Long - 21.00 each
 Cat. No. 875-014 Kopp Right Long - 21.00 each
 Cat. No. 875-015 Kopp 45° Right - 21.00 each
 Cat. No. 875-016 Kopp 45° Left - 21.00 each

SPECIAL

1 set 6 Kopp
Cat. No. 875-010
One of each style
 List \$160 for \$150

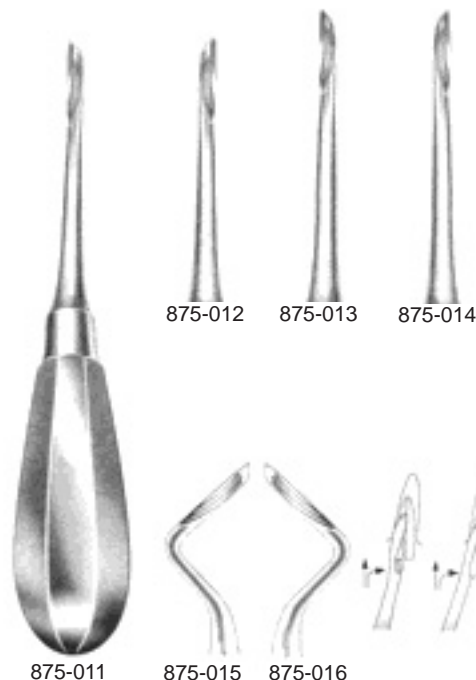

Claw Elevators

- Razor sharp tip for piercing mandibular molar interseptal bone and lifting root tip.
- Each instrument has a right and left side to work both mesially and distally.
- Bent or dull tip instruments should be replaced.

Cat. No. 875-021 Right Hand Claw Elevator - 35.00 each
 Cat. No. 875-022 Left Hand Claw Elevator - 35.00 each

SPECIAL

1 set 2 Claw
Cat. No. 875-020
One of each style
 List \$70 for \$63



Bone Mill

- A compact quality bone mill that is easy and safe to operate and ECONOMICAL.
- The all stainless steel parts are simple to assemble and completely disassemble in order to autoclave.
- Reduces cortical and cancellous bone to desired particle size by varying gentle hand grinding pressure.
- Resulting bone particles are used for supporting implants, augmentation, and repair.

Cat. No. 875-040 - 350.00 each

SPECIAL

Bone Mill
Totally
Autoclavable

